

# U0007

## FUNCTIONS



Backrest adjustment : 0 - 68°



Thigh rest adjustment : 0 - 27°



Foot rest mechanical adjustment : 0 - 22°



Vascular position



Fowler position



Cardiac position



Trendelenburg position : 0 - 16°



Anti - trendelenburg position : 0 - 16°



Mattress platform minimum height : 49 cm (150 mm with castor)



Mattress platform maximum height : 79 cm (150 mm with castor)


# HOSPITAL BED 4 MOTORS WITH SCISSOR

## OPTIONS

### GROUP 1 / CASTOR OPTIONS



### GROUP 2 / SIDE RAIL OPTIONS

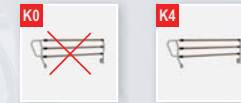


### GROUP 3 / ELECTRICAL ACCESSORY OPTIONS



### GROUP 4 / BACKREST MOTOR OPTIONS




### GROUP 5 / LEGREST MOTOR OPTIONS

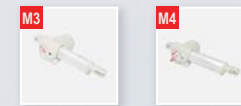

### GROUP 6 / MATTRESS PLATFORM OPTIONS



### GROUP 7 / OTHER ACCESSORY OPTIONS

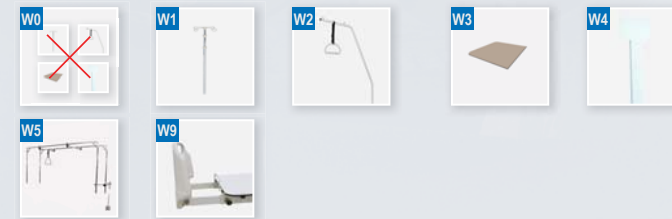

## DESCRIPTION OF OPTION CODES

<b>T16</b> CASTOR TWO WITH BRAKE (125 MM)	<b>M3</b> LEGREST MOTOR (4000 N)
<b>T14</b> CASTOR CENTRAL & DIRECTIONAL LOCK (125 MM)	<b>M4</b> LEGREST CPR MANUAL CONTROL MOTOR (2000 N)
<b>T19</b> CASTOR CENTRAL & DIRECTIONAL LOCK (150 MM)	<b>Y1</b> REMOVABLE 4 PARTED MESH PLATFORM.
<b>K0</b> WITHOUT SIDE RAIL	<b>Y3</b> REMOVABLE DOUBLE WALLED PE BED PLATFORM
<b>K4</b> WITH ALUMINIUM SIDE RAIL	<b>W0</b> WITHOUT OTHER ACCESSORIES
<b>E0</b> WITHOUT ELECTRICAL ACCESSORIES	<b>W1</b> IV STAND : ONE HEIGHT ADJUSTMENT SYSTEM
<b>E1</b> BATTERY RECHARGEABLE LEAD GEL TYPE. OUTPUT VOLTAGE: 24 V DC. CAPACITY: 1,2 Ah.	<b>W2</b> LIFTING POLE : BEARING CAPACITY 75 Kg
<b>E2</b> NURSE CONTROL PANEL. PROTECTION CATEGORY : IP 66. 15 FUNCTION KEY. UNDER BED LIGHTING SYSTEM	<b>W3</b> X-RAY AVAILABLE FOR BACK SIDE WITH CASSETTE HOLDER AND COMPACT LAMINATE
<b>E3</b>	<b>W4</b> TELEVISION STAND FOR 17" LCD TV
<b>M1</b> BACKREST MOTOR (4000 N)	<b>W5</b> TRACTION APPARATUS
<b>M2</b> BACKREST CPR MANUAL CONTROL MOTOR (2000 N)	<b>W9</b> BED EXTENSION : 20 cm.

Must select only one of them.

Do not select the others if you select without accessory. If you select the others you can select all or 1-2-3 ... of them.

## GENERAL SPECIFICATIONS

- Easy to clean four-piece bed platform.
- Height, trendelenburg, reverse trendelenburg, head and foot movements with the help of an electrical motor.
- Ergonomic, durable and easily cleanable aluminum rails system with maximum safety designed in accordance with the standards of EN 60601-2-52.
- Head and foot platforms with automatic sliding system.
- Easily removable and installable, hygienic, high resistant head board system designed of compact laminate and aluminum.
- Shock resistance plastic bumper wheels.
- Metal components dyed with anti-corrosion electrostatic powder paint after being treated with phosphate.
- Profile frame with easy removable PE plastic covers.
- Optional laminate color choice.<sup>(RLX)</sup>

## TECHNICAL SPECIFICATIONS

■ External dimension	: 102 x 222,5 cm
■ Width outside side rails	: 98 cm
■ Mattres platform	: 89 x 195 cm
■ Bed extension	: +27 cm
■ Recommended mattress dimension	: 85 x 194 x 12 cm
■ Autoregression (Backrest / Legrest)	: 8 cm / 4 cm
■ Safe Operating Load (Inc. Mattress and Accessories)	: 230 kg
■ Weight	: ~ 127 kg
■ Volume	: 1,23 m <sup>3</sup>

## ELECTRICAL SPECIFICATIONS

■ Power Input	: 230 V AC 50 Hz Max 1 A
■ Power Output	: 24 V DC
■ Electrical Current Protection	: Class 1, Type B
■ Liquid and Dust Protection	: IP54 – IP66



\* RL : Laminate color options page 82.