

FUNCTIONS



Backrest adjustment : 0 - 65°



Thigh rest adjustment : 0 - 48°



Foot rest mechanical adjustment : 0 - 25°



Vascular position



Fowler position



Trendelenburg position : 0 - 12°

OPTIONS

GROUP 1 / CASTOR OPTIONS



125 Kg





125 Kg

GROUP 2 / SIDE RAIL OPTIONS



GROUP 3 / OTHER ACCESSORY OPTIONS

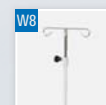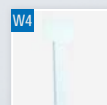

DESCRIPTION OF OPTION CODES

T16 CASTOR TWO WITH BRAKE (125 MM)

T4 METAL CASTOR TWO WITH BRAKE (125 MM)

K0 WITHOUT SIDE RAIL

K1 WITH METAL SIDE RAIL

K3 WITH METAL SIDE RAIL

Must select only one of them.

W0 WITHOUT OTHER ACCESSORIES

W2 LIFTING POLE : BEARING CAPACITY 75 Kg

W4 TELEVISION STAND FOR 17" LCD TV

W5 TRACTION APPARATUS

W8 IV STAND : ONE HEIGHT ADJUSTMENT SYSTEM

Do not select the others if you select without accessory. If you select the others you can select all or 1-2-3 ... of them.

GENERAL SPECIFICATIONS

- Easily cleanable metal bed platform designed for the ventilation of the bed.
- Head and foot movements with the help of an electrical motor.
- Trendelenburg movement by the help of shock absorber mechanism at foot side.
- Ergonomic, durable and easily cleanable metal rails system with maximum safety.
- Melamine coated chipboard head and foot panels.
- Removable and installable head system that can be also used as a resulation board if necessary.
- Shock resistance plastic bumper wheels.
- Castor system with two-castor brake in a diameter of 125 mm.
- Metal components dyed with anti-corrosion electrostatic powder paint after being treated with phosphate.
- Optional melamine coated chipboard color choice.^(RSX)

TECHNICAL SPECIFICATIONS

■ External dimension	: 105 x 204 x 50 cm
■ Recommended mattress dimension	: 85 x 190 x 12 cm
■ Safe Operating Load (Inc. Mattress and Accessories)	: 230 kg
■ Weight	: ~ 86 kg
■ Volume	: 1,17 m ³

ELECTRICAL SPECIFICATIONS

■ Power Input	: 230 V AC 50 / 60 Hz 50W
■ Power Output	: 24 V DC



* RS : Melamine color options page 83.